




- PBM (Photo Bio-Modulation) Applications
- PDT (Photo Dynamic Therapy) Applications.
- Portable
- Battery Operated


Benefits

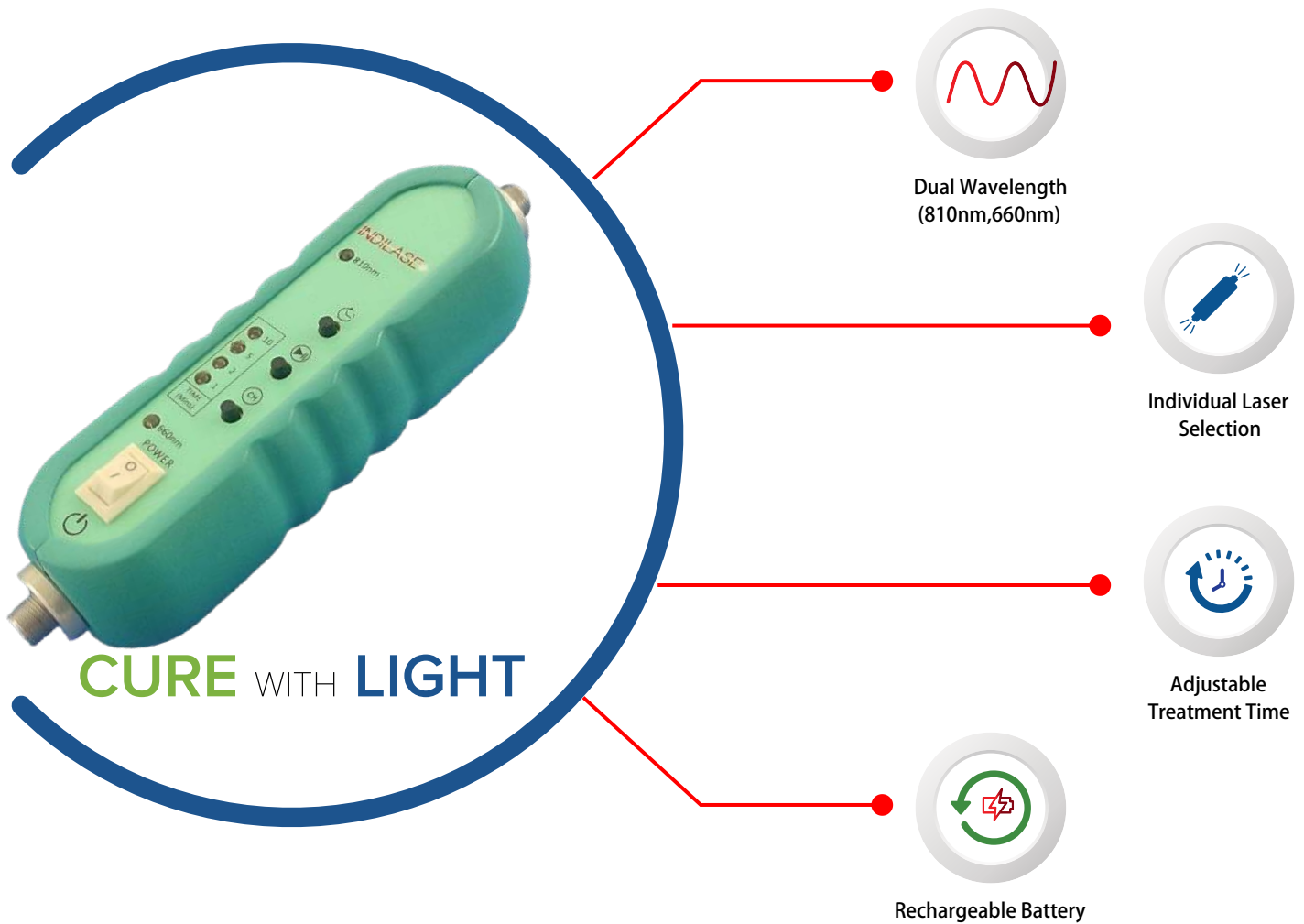

Pain Free



Quick Healing

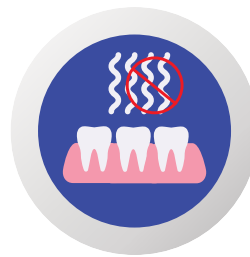

Non Thermal

Soft and Hard Tissue Applications of Low Level Laser Therapy (LLLT)

Soft Tissue Applications	Wavelength
Herpetic Lesions	660 nm / 810 nm
Aphthous / Traumatic Ulcers	660 nm / 810 nm
Post Oncology Mucositis	660 nm / 810 nm
Post Extraction Socket / Post Trauma Sites	660 nm / 810 nm
Trigeminal neuralgia	810 nm
Periodontal Pocket Disinfection / Periodontitis	660 nm (Use with Methylene Blue Dye) / 810 nm (Use with Toluidine Green Dye)
Edema	810 nm
Sinusitis	810 nm
Gag Reflex / Nausea	
Post Surgical Pain	660 nm / 810 nm
Mouth Opening	660 nm / 810 nm

Product Specifications

Wavelength and Power	660 nm / 100 mW 810 nm / 200 mW
Operation Mode	Continuous Output
Treatment Time (Minutes)	1, 2, 5, 10.
Dimensions (L x W x H)	(140 x 25 x 40)mm
Weight	110 Gram
Warranty	1 Year

Hard Tissue Applications	Wavelength
Temporomandibular Disorders	810 nm
Dental Hypersensitivity	810 nm
Stimulatory effect on root development	810 nm
Implants	810 nm

Accessories



Bio-Stimulation Attachment 8mm



TMJ Attachment 20mm



Safety Goggle



Eye Protector



Carry Case

MEDSOL

#304A / G2, Brindavan Nagar, T.N.H.B, Phase-VII,
Hosur-635109, Tamilnadu, India.

Customer Care : +91 63829 67677 / 04344 247499

Email : info@indilase.co.in

www.indilase.com

www.indilase.com